
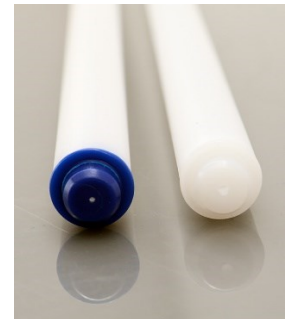


1. Part Number: 8051H-450S

2. Photograph:

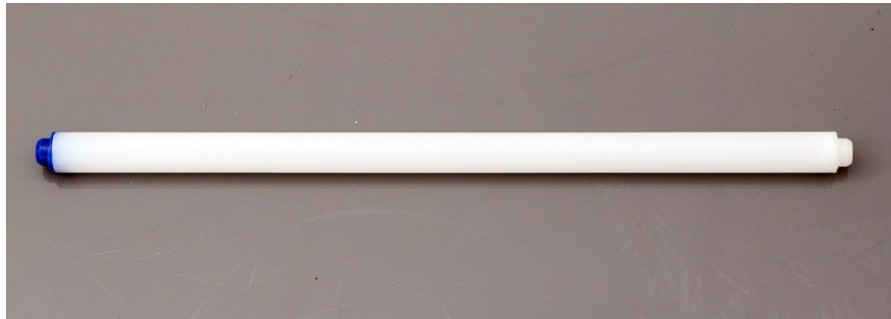

3. Product Information:

Description: SteriWare® Jumbo Pipette 450mm - Sterile

Use: Single use, disposable liquid sampler

Material of Construction: High Density Polyethylene (HDPE) - virgin

- Conforms to FDA 21 CFR 177.1520
- Conforms to EU Regulations 10/2011
- Conforms to EC Regulations 1935/2004
- BSE/TSE free

A blue dye is added to the Tip. The blue dye certified as being FDA acceptable and does not contain any materials of animal origin.

Method of Construction: Ultrasonic welding

Nominal Length: 420mm

Diameter of Outer Tube: 21mm

Diameter of hole in tip: Handle End (Blue) = 1.5mm
Product End (Natural) = 3mm

Nominal Weight of Sampler: 50g

Moulding Environment: Class 100,000 Cleanroom

Assembly & Primary Packing Environment: Class 10,000 Cleanroom

Individually Bagged? Yes (heat sealed PE/PET bag)

Method of Sterilisation: Gamma Irradiation (25 kGy minimum dose)

Number of Samplers per Box: 50

Recommended Storage Conditions: Dry and ambient temperature

Shelf Life: 5 years from the date of manufacture

ATEX: This product has not been ATEX tested

4. Change Control:

Sampling Systems operates a policy of continual improvement.

To be kept informed of all major changes to any product in the SteriWare range please email your contact details to change.control@sampling.com

5. Country of Origin:

Made and packed in the UK.

6. Repeated Use:

This is a single-use, disposable sampler. Repeated use will lead to failure.

7. Approval:

Signed:



Name: Kuldeep Flora

Position: Quality Manager

Date: 11-Jul-24