


**1. Part Number:** A416-20

**2. Photograph:**

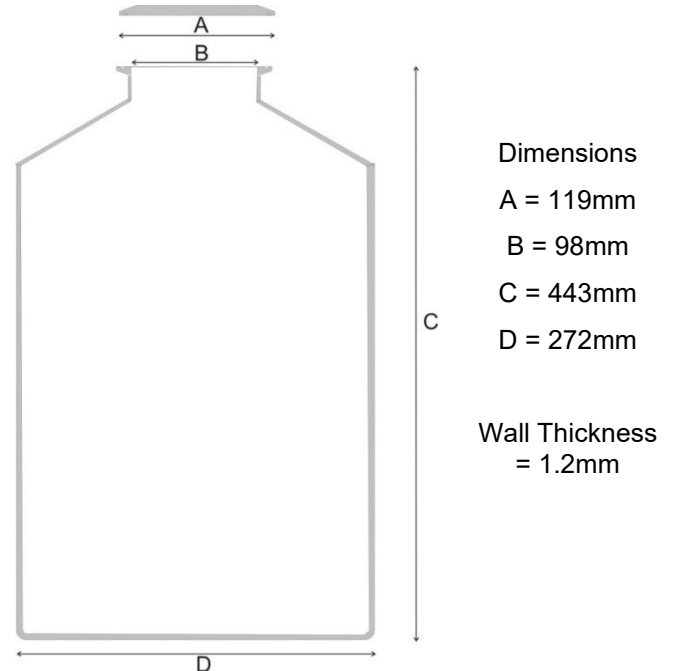

### 3. Product Information:

Description:	Stainless Container 20L
Materials of Construction:	<ul style="list-style-type: none"> <li>- Body: 316L stainless steel</li> <li>- Lid: 316L stainless steel</li> <li>- Clamp: 304 stainless steel</li> <li>- Gasket: Silicone (FDA acceptable grade)</li> </ul>
Method of Construction:	All welds ground and polished Crevice free interior
Surface Finish:	Better than 0.5 microns Ra
Nominal Volume	20L
BSE/TSE statement:	All polishing compounds used in the manufacture of this product are of vegetable origin.
Weight:	5.1 Kg (body + lid + gasket + clamp)
Recommended Storage Conditions:	Dry and ambient temperature
Pressure Rating:	None – these containers are not pressure rated vessels
Notes:	<ul style="list-style-type: none"> <li>- The neck of the container is fitted with a 4" diameter ferrule which conforms to BS4825-3</li> <li>- These containers are primarily designed for storing powders</li> <li>- They are not suitable for the storage or transportation of beer</li> </ul>

## 4. Spare Gaskets:

Supplied with a 4" silicone gasket as standard. Gaskets are consumable parts and will need to be inspected and changed at the first signs or wear.

Part No.	Description
K201-4	Gasket 4" silicone – FDA acceptable grade
K202-4	Gasket 4" EPDM – FDA acceptable grade
K203-4	Gasket 4" Viton – FDA acceptable grade

## 5. Approval:

Signed:



Name: Kuldeep Flora

Position: Quality Manager

Date: 14-Apr-26